ASSEMBLY INSTRUCTIONS


Pedal Exerciser with electronic display

- NL** *Fietstrainer met elektronisch display*
- FR** *Pédalier avec afficheur numérique*
- DE** *Fahrradtrainer mit elektronischem Display*
- ES** *Pedalier con pantalla electrónica*
- PT** *Aparelho de pedalar com mostrador eletrónico*

Aids to daily living

Able2

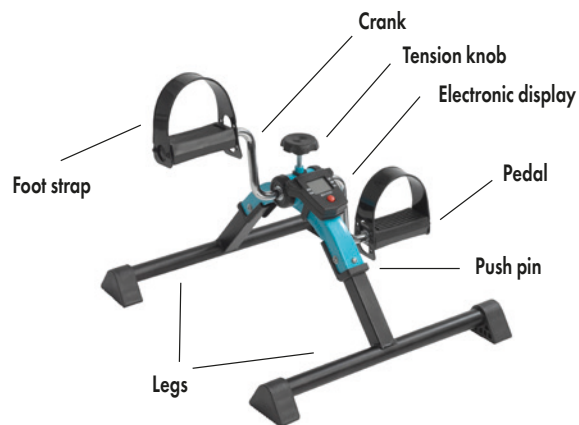1. Remove contents from box.
2. Pull legs in opposite direction until locking button "snaps" into locking receiver. (see figure A)
3. Attach electronic display to the pedal exerciser by sliding the loop over the tension knob around the frame (see figure B)

Attention: Only attach the multi function display to this adhesive label!

4. Insert tension knob on top of pedal exerciser and turn clockwise to tighten. The tension knob will increase/decrease the amount of tension on the pedal exerciser. Tighten/loosen to a level that feels comfortable.
5. You can set the width of the foot strap through the 2 adjustment holes. (see figure C)

figure A

figure B

figure C

Attention: Do not stand on the pedal exerciser! Do not put too much weight on the pedals, otherwise the pedal exerciser may be damaged.

FOLDING

To fold, depress push pin and fold leg inward. Repeat for second leg. (see figure D).

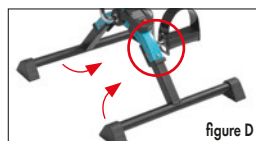
figure D

folded

FUNCTION AND DISPLAY DESCRIPTION

1. **MODE (Red Button):** this control allows you to select a particular function. To Start, press and hold the Red Button for 4 seconds. The LED will display the different functions. If you do not select a function, the unit will automatically shut off.
2. To select a function, wait until the 4 seconds set up is complete. Push the Red Button and the arrow will move through the functions. When the desired function is selected, simply use the pedal exerciser and the unit will analyse the exercise.
3. **COUNT (CNT):** counts the number of steps.
4. **TIME:** indicates duration of training.
5. **CALORIE (CAL):** indicates the amount of calories burned.
6. **REPS PER MINUTE (RPM):** automatically calculates the number of repetitions per minute.
7. **SCAN:** automatically reviews the results of all functions.

To reset the display to zero, press the Red Button (RESET) 3 seconds.

BATTERY REPLACEMENT

If the display does not come on or reads incorrectly, the battery requires replacement. To replace the battery slide the lid on the back of the display down, replace the underlying battery and set the display back on.

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Valentine International Ltd. Taiwan
Able2 Europe B.V.
Able2 UK Ltd.